


# ENFit Enteral System Connector Changes


The new design standard impacts the entire enteral feeding system

NUTRITION END

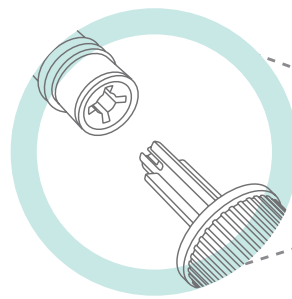


**CONNECTOR (FINAL)**  
[In place since 2012]


PATIENT-ACCESS END

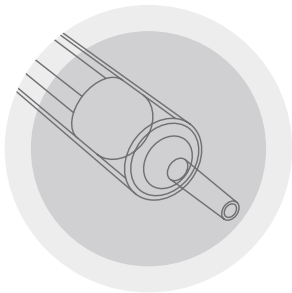

**SYRINGE (CURRENT)**

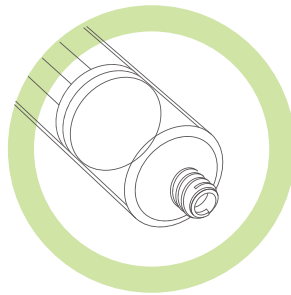

OR

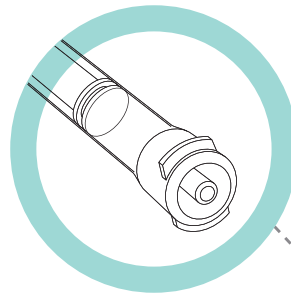


**SYRINGE (Standard Tip)**


Syringes to administer medicine, flush, hydrate, or bolus feed through enteral tubes will now require a precise enteral-specific fitment.

**SYRINGE (Low Dose Tip)**

To ensure small volume dosing accuracy, syringe sizes of 5mL or smaller may require an ENFit Low Dose Tip.

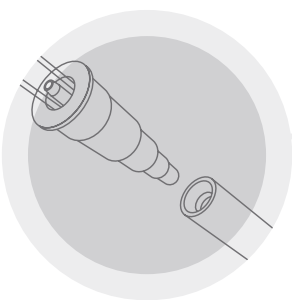


**FEEDING TUBE (CURRENT)**

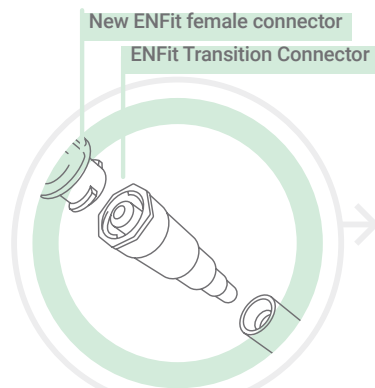

**TRANSITION SET (TEMPORARY)**

Allows fitment to current feeding port until new ENFit enteral feeding tubes are available.

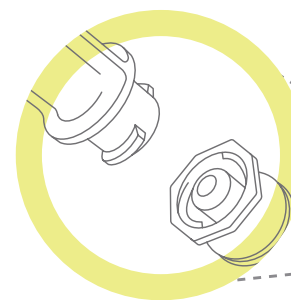

**FEEDING TUBE (FINAL)**

Changing from male—the stepped or Christmas tree connector—to the new ENFit female connector. The feeding tube port for the administration set will change from female to male.

**Note:** Speak to your supplier representatives for availability, timing and indications for use of ENFit administration sets, syringes and feeding tubes.



For more information and to sign up for email updates, visit [StayConnected.org](http://StayConnected.org).

v.2 ©2016 GEDSA